

# Bowie & Medical

BODY CARE SYSTEMS

## Bowie IPL 2.0



IPL equipment  
with PPT Technology

Bowie Medical B.V. Parallelweg 7d 5349 AD Oss The Netherlands  
Tel.: 0031 (0)412 69 01 16 Fax 0031(0)412 69 01 17  
Info: [www.bowiemedical.nl](http://www.bowiemedical.nl) E-mail: [info@bowiemedical.nl](mailto:info@bowiemedical.nl)


Bowie-IPL 2.0: design and engineering by Bowie Medical BV - The Netherlands  
Bowie Medical BV in The Netherlands is the worldwide exclusive supplier of Bowie IPL equipment.

## Bowie IPL

The Bowie IPL 2.0 is a powerful water-cooled system with a working capacity of up to 50 Joules per cm<sup>2</sup>. The system is equipped with many features, such as variable pulse-dosing, adjustable pulse-sequences, rapid pulse-repetition. The system operates with extremely low costs per pulse. It has full Medical CE Certification. The Bowie IPL can be used for permanent hair removal for all skin types and hair types. The device can also be used to treat: Rosacea, Acne, Pigmentation, Skin Blemishes and Skin-Rejuvenation.

## PPT Technology

The new BOWIE IPL 2.0 uses the latest "Perfect Pulse Technology" (PPT). The PPT technology provides a controlled perfect pulse-shape (block-shape) with constant and equal delivery of power during the full pulselength. PPT provides better reproducible results as well as gentler, more comfortable treatments. See chart below:

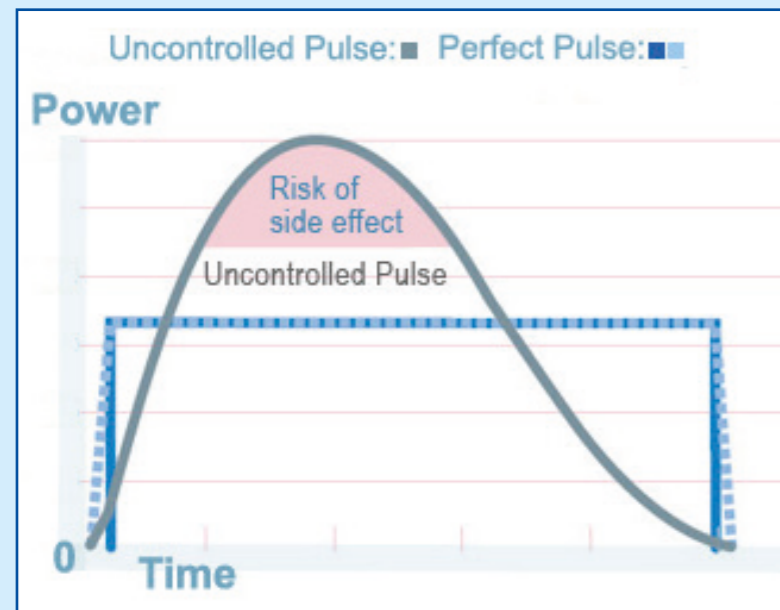

## Many Advantages

The new Bowie IPL 2.0 is equipped with a new intuitive software system with a full color touch screen. The system works with easily replaceable filters. The Bowie IPL can be used for permanent hair removal for all skin types and hair types. The device can also be used to treat: Rosacea, Acne, Pigmentation, Skin Blemishes and Skin-Rejuvenation. The powerful Bowie IPL 2.0 has full Medical CE Certification.



## Specifications

- With European Medical CE Approval (CE Certificate)
- Full color intuitive touchscreen control panel
- Menu in different languages
- Very high working (medical) capacity (up to 50 J / cm<sup>2</sup>)
- Spot Size: 15mm x 50mm (7.5 cm<sup>2</sup>)
- Easily replaceable filters
- High Performance water cooling
- Adjustable skin cooling
- Very fast pulse repetition (0,2 - 2,0 Hz)
  - Variable pulse setting
  - Pulse Series Programming (1-99 pulses)
- Ergonomic handpiece
- Both foot & hand control
- With latest PPT technology inside



### Available filters for Bowie IPL:

430nm - 530nm - 590nm - 640nm - 690nm - 750nm  
(standard 2 filters included with purchase: 530nm and 640 nm)

**Dimensions:** 600x350x430 mm - Weight: 28 kg

Peak Power: 1200 Watt - Pulse Duration: 0,5 - 99 ms

Pulse Delay: 1-99 ms

## Service & Training

Bowie Medical BV will provide full technical and educational (training) services.

## About Bowie Medical

BOWIE MEDICAL BV in THE NETHERLANDS  
Has developed the new Bowie IPL 2.0 System for permanent hair-removal and other skin therapies. The system is the result of excellent Dutch design and engineering.

